

TUTORIAL

# Setting up Prisma Touch Screens

Prisma 20A/20C

LÖWENSTEIN  
medical

 **santé**  
NEXT-GENERATION SOURCING



# CONTENTS

Clinical Menu	1
Clinical Info Menu	6
Humidifier	7
RAMP (softSTART)	8
Notes	9

# CLINICAL MENU

1. Use a firm touch and hold when navigating the Prisma touch screen devices, the touch screen is thick and durable and therefore not as sensitive as a smart phone.




2. Open the clinical menu:
  - Hold down the info button for > 4 seconds ----> 
  - A confirmation message appears, press "Login"

3. Exit the clinical menu:
  - Press the Logout button ----> 
  - A confirmation message appears, press "Logout"

4. Clinical menu home/opening screen (see right)



1. Therapy mode rapid access
2. Main menu button
3. Ventilation parameters rapid access
4. Clinical info menu (explained later in this document)

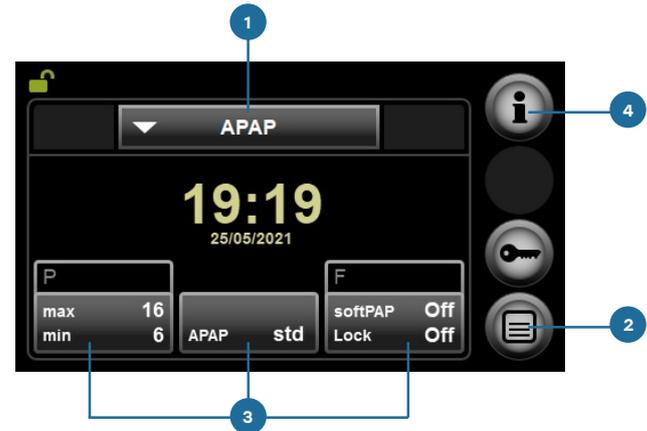

## Tutorial


Setting up Prisma 20A/20C Touchscreens

- Per “3”, you can immediately select the parameter you wish to change and adjust by either sliding the scale, pressing the desired setpoint on the scale, and/or using the “+” and “-” buttons. You can change as many parameters as available on the page (see below) and once satisfied with all adjustments press the green tick in the bottom right to confirm.

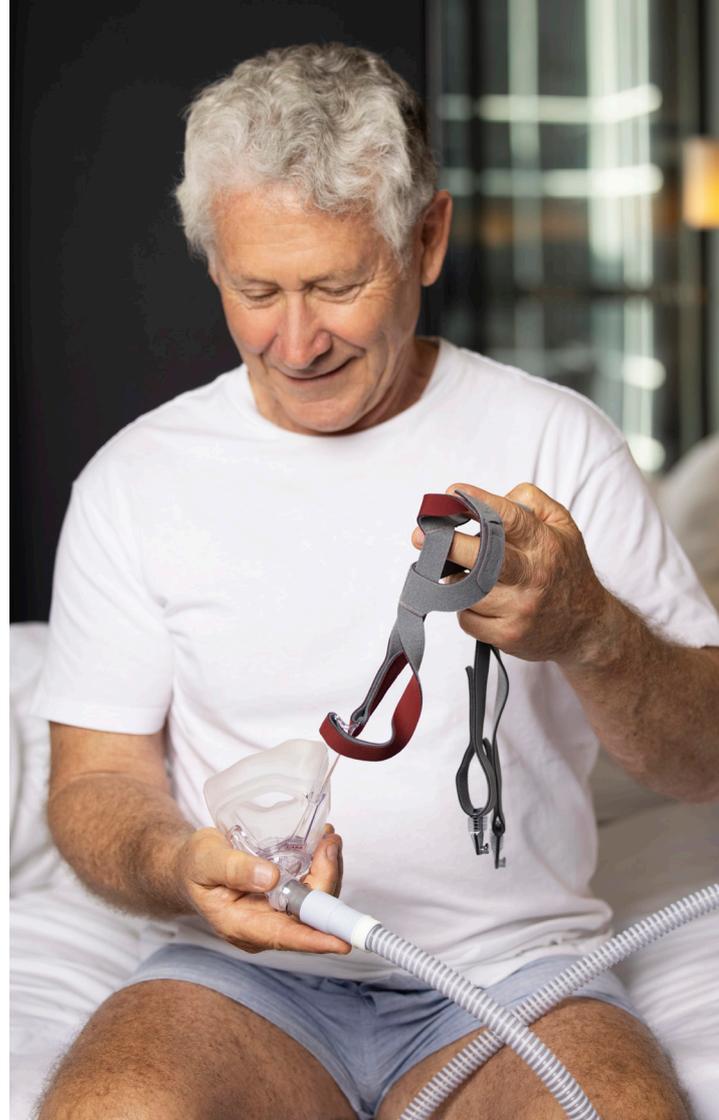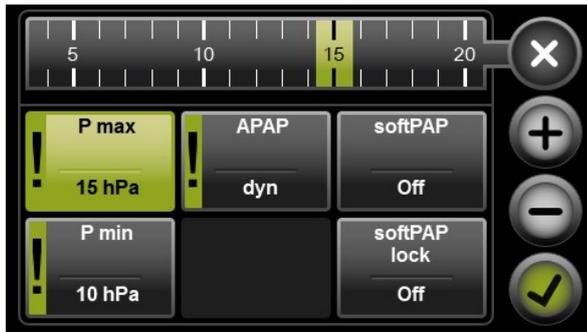


6. Navigate to main menu screen by pressing the button in the bottom right (“2. Main Menu Button”).

The following options appear.

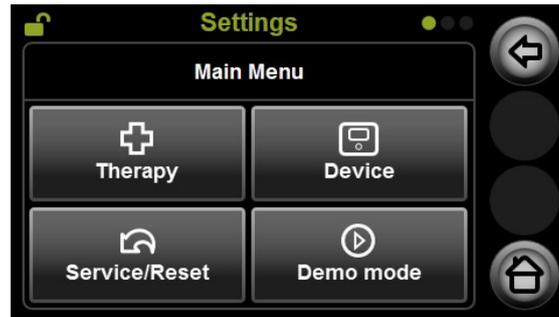

When **Therapy** is pressed there are 4 sub sections:

1. **Ventilation parameters** (options per rapid access in point 3.)
2. Alerts (any required alarms)
3. Comfort (set AutoStart and/or AutoStop, and ramp settings [called softSTART])
4. Accessories (option to adjust tube diameter and specify if a bacteria filter is being used)

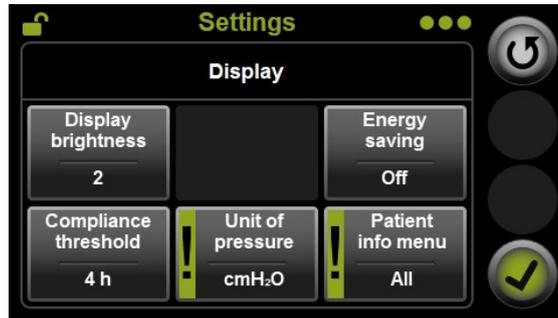

When **Device** is pressed there are 5 sub sections:

1. **Time** (set date & time, in clinical menu the hour setting can also be changed if required)
2. **Signal Data** (no. of days of signal data, recommend leaving this at 14 days)
3. **Language** (ability to change language for patient vs. clinical menu)
4. **Display** (ability to change various display settings, recommend changing “hPa” to “cmH20” and “Patient Info Menu” to “All” [this will allow patient to see usage time, median leakage, and AHI per hour in the patient menu]). Once selected, use the green tick to confirm – see below image.
5. **Audio** (you can remove/change key tone volume as well as change alarm volume)

When **Service/Reset** is pressed there are 3 options to then press:

1. **Reset to factory settings** (generally unused)
2. **Delete therapy data** (i.e. to clear data when transferring device to a new patient)
3. **Maintenance counter** – e.g. for filter replacement

**Demo mode** is generally unused, this runs through a cycle of screen shots.



## Tutorial

Setting up Prisma 20A/20C Touchscreens

# CLINICAL INFO MENU

Press “i” in top right once whilst in clinical menu to open clinical info menu. Further information can be gathered by pressing 1, 2, or 3 as labelled below.




- 1 Selection  
Therapy period
- 2 Therapy quality  
Detailed information, e.g.:
  - obstructive / central AHI
  - RERA
  - Snore
  - Periodic breathing
- 3 Leakage  
Detailed information



# HUMIDIFIER

When therapy is running, if the humidifier is on the symbol will light up green with a number (as per image on left). Use “+” or “-” to change humidity level. Press the humidifier symbol once to turn off and it will go grey (image on the right). Press again to reactivate.

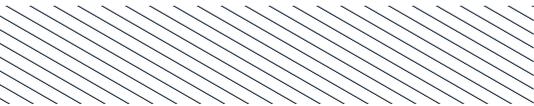


## Tutorial

Setting up Prisma 20A/20C Touchscreens

# RAMP (softSTART)

Similar to humidifier, ramp is on when the symbol is lit up green with a number which is the remaining ramp time (as per image on the left). Press the ramp symbol once to turn off and it will go grey (as per image on the right). Press again to reactivate.

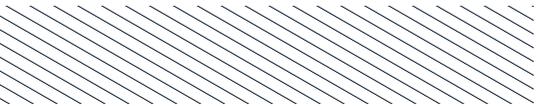



# LÖWENSTEIN medical



[lowensteinmedical.com.au](http://lowensteinmedical.com.au)